Eta Body Balance

osobní váha

ETA777990000

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Dear customers,
Thank you for your purchase and interest into our personal scale. Use this manual together with the maual for the scale.This manual is focused only on the application and its options. It is not focused on the usage of the scale.

Illustrations and translations of texts or of activating elements applied in the application may slightly differ from the illustrations and descriptions specified in this manual. Therefore, don't get scared if there are differences in your application compared with the illustrations in the manual. The illustrations provided may also cover different measurements, so don't be afraid if you see a discrepancy between the results.

CONTENT	
BASIC STRUCTURE OF THE MAIN HOME SCREEN OF THE APPLICATION	24
DEVICE AND USER MANAGEMENT	25
SCALE SETUP AND MANAGEMENT	26
USER MANUAL FOR THE Eta Body Balance	26
APPLICATION	27
1. APPLICATION INSTALLATION	27
2. ACCOUNT REGISTRATION AND LOGIN	27
3. DEVICE PAIRING	28
4. BASIC ADVICE AND TIPS FOR USE	29
5. WEIGHING AND ANALYSIS	29
6. MEASUREMENT RESULTS	30
7. SHARING RECORDS AND MEASUREMENT RESULTS	31
8. THIRD-PARTY APPLICATIONS	32
9. QUESTIONS AND ANSWERS	32
10. FEEDBACK	32

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

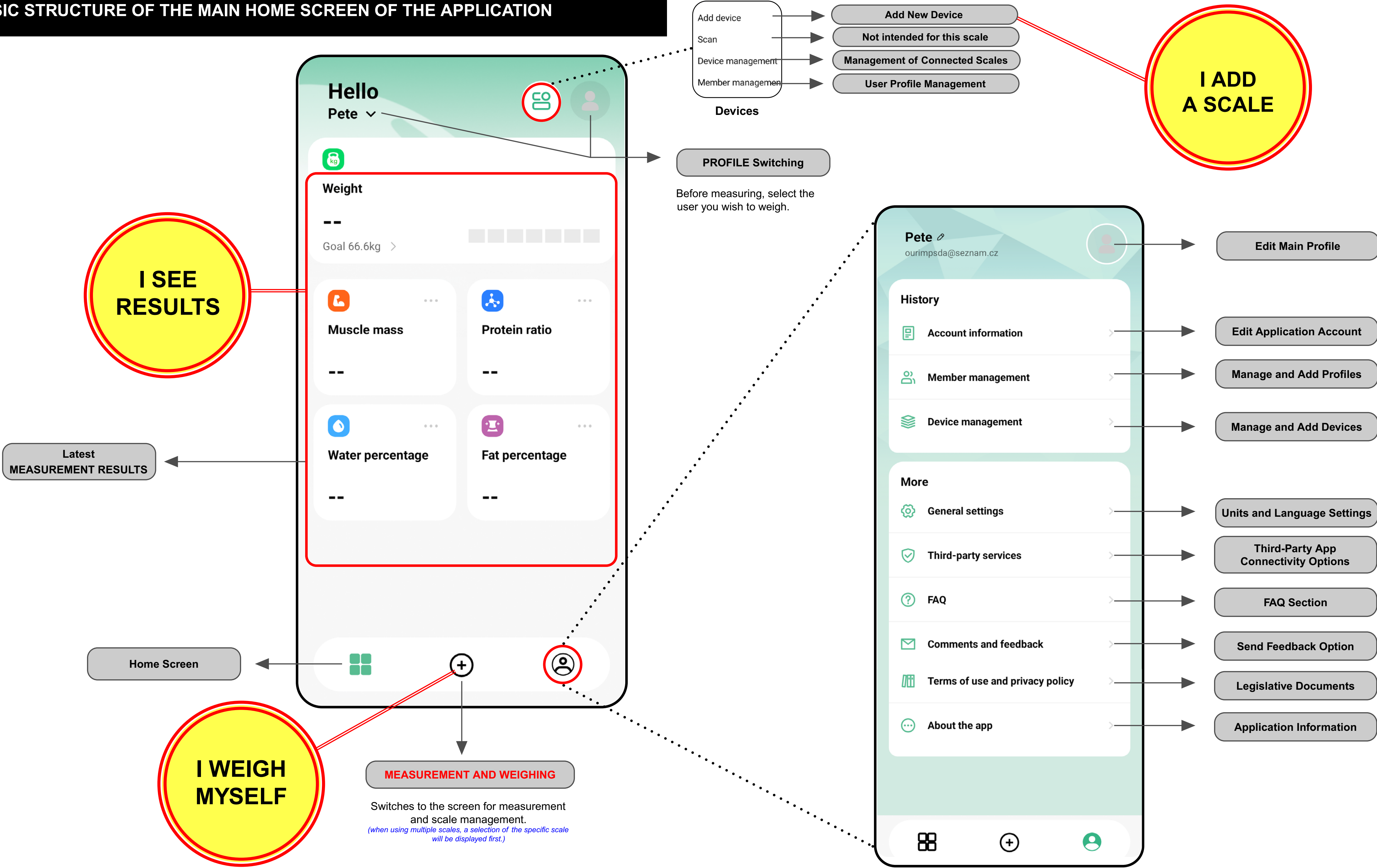
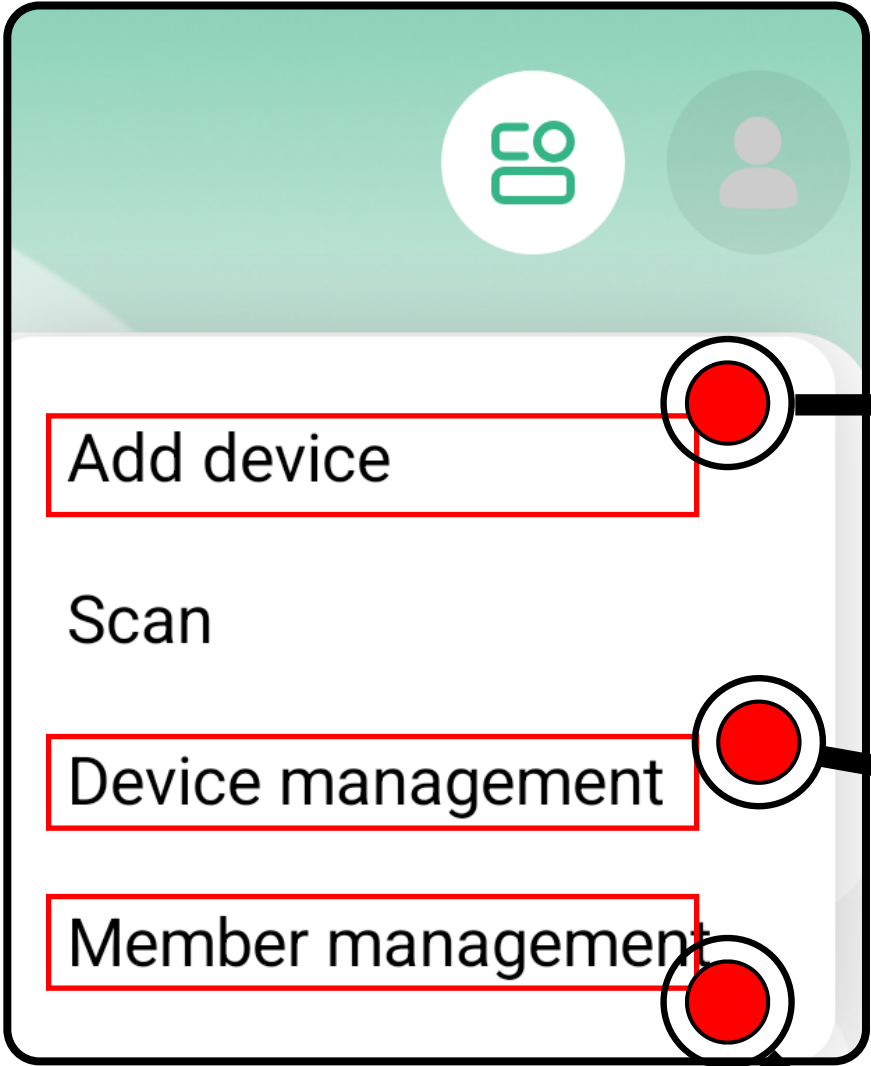
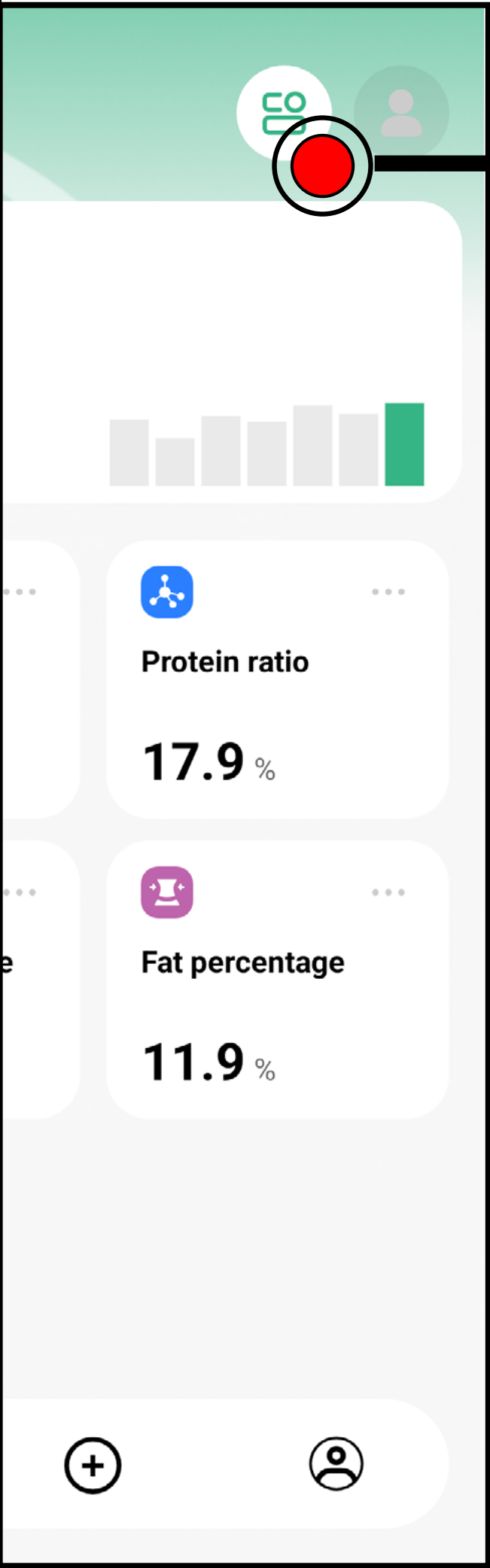
ADD
A SCALE

Selecting the type of scale and subsequent **PAIRING**

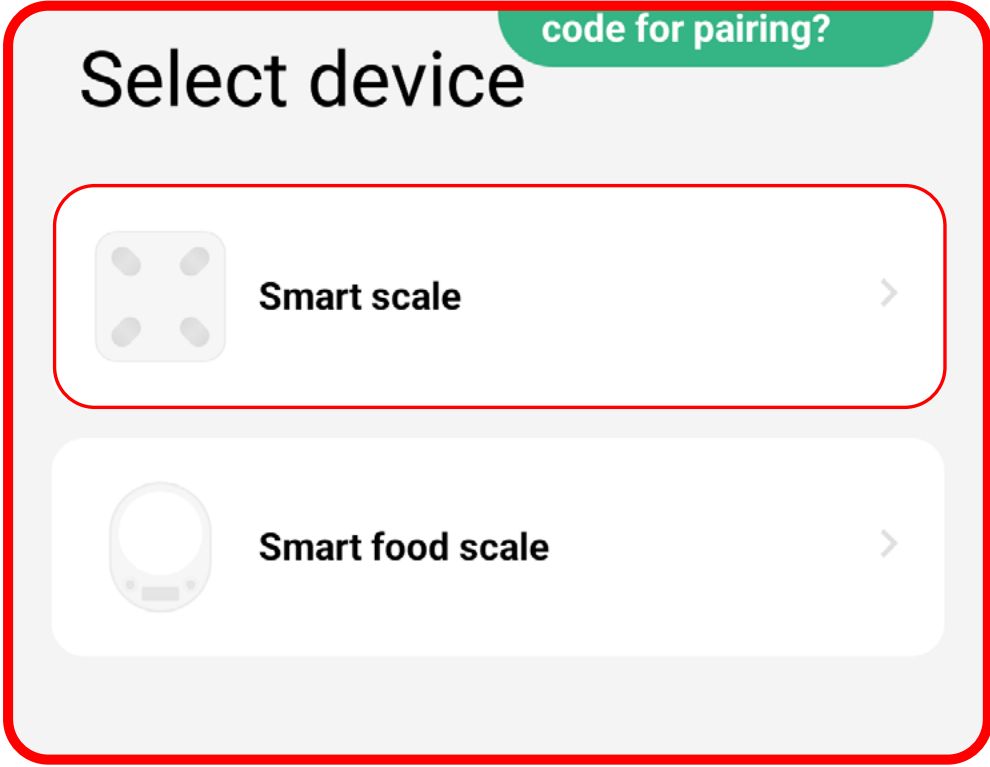
CREATE A
PROFILE

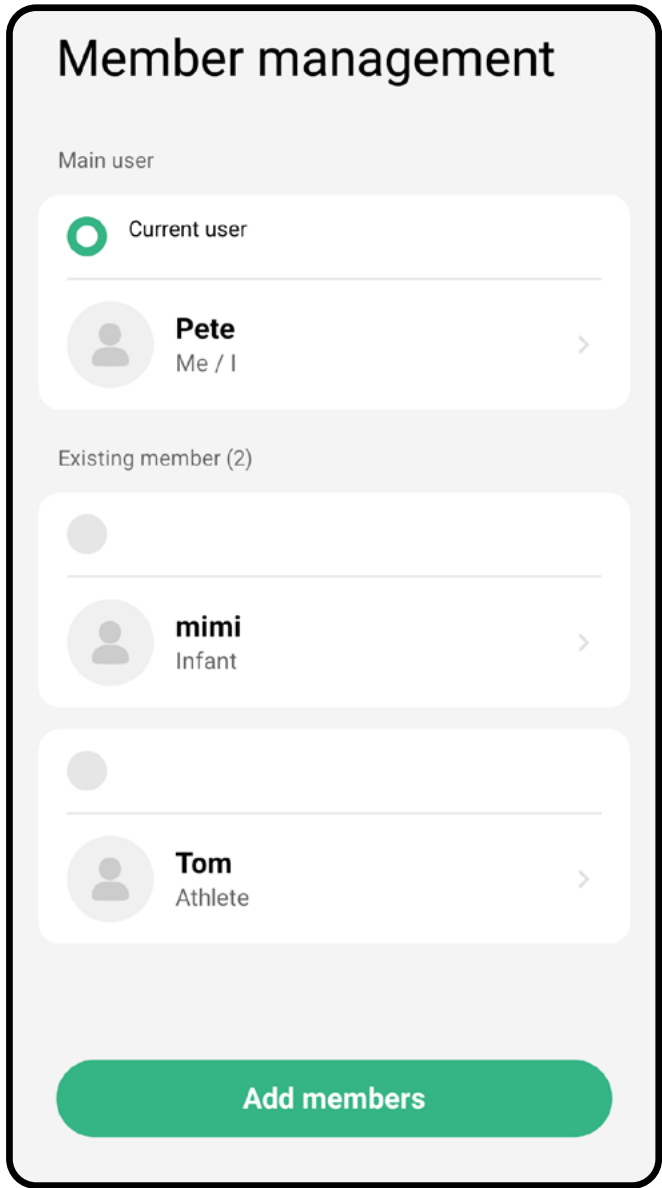
Managing and adding **USER PROFILES**
with the option to edit them.

(To delete a profile, press and hold it — a delete option will then appear.)

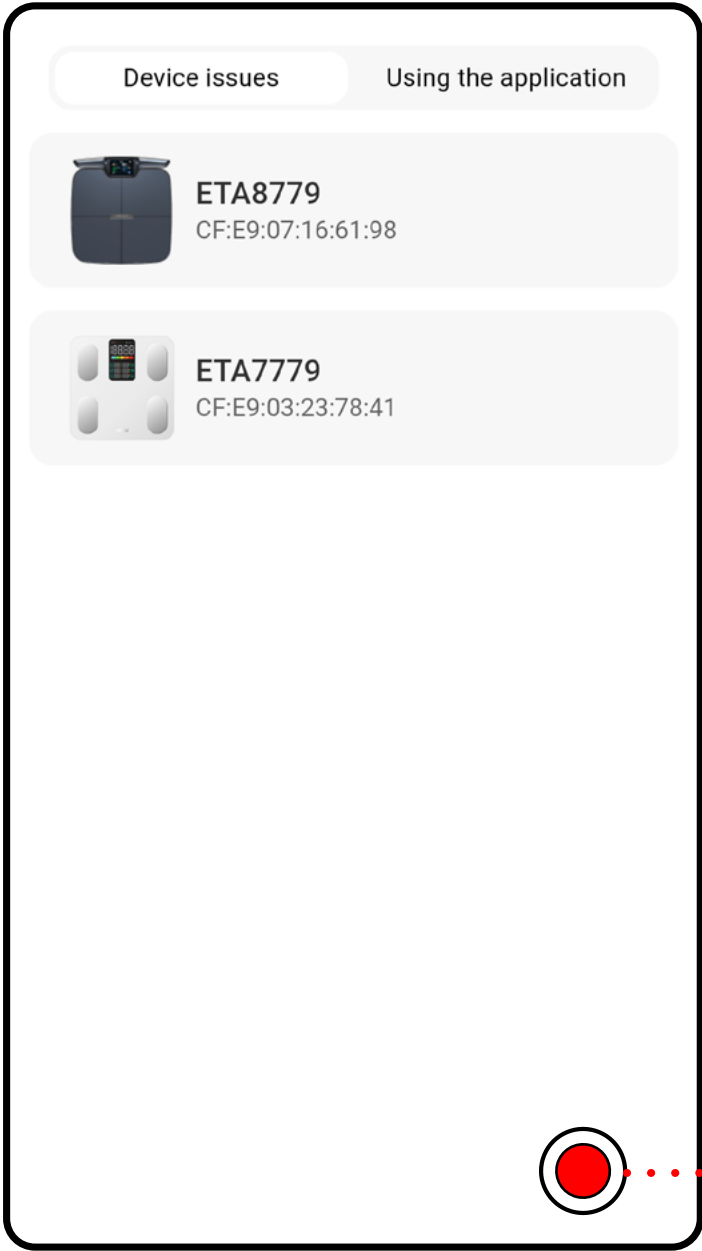Management of all paired devices,
with the option to add a new device.

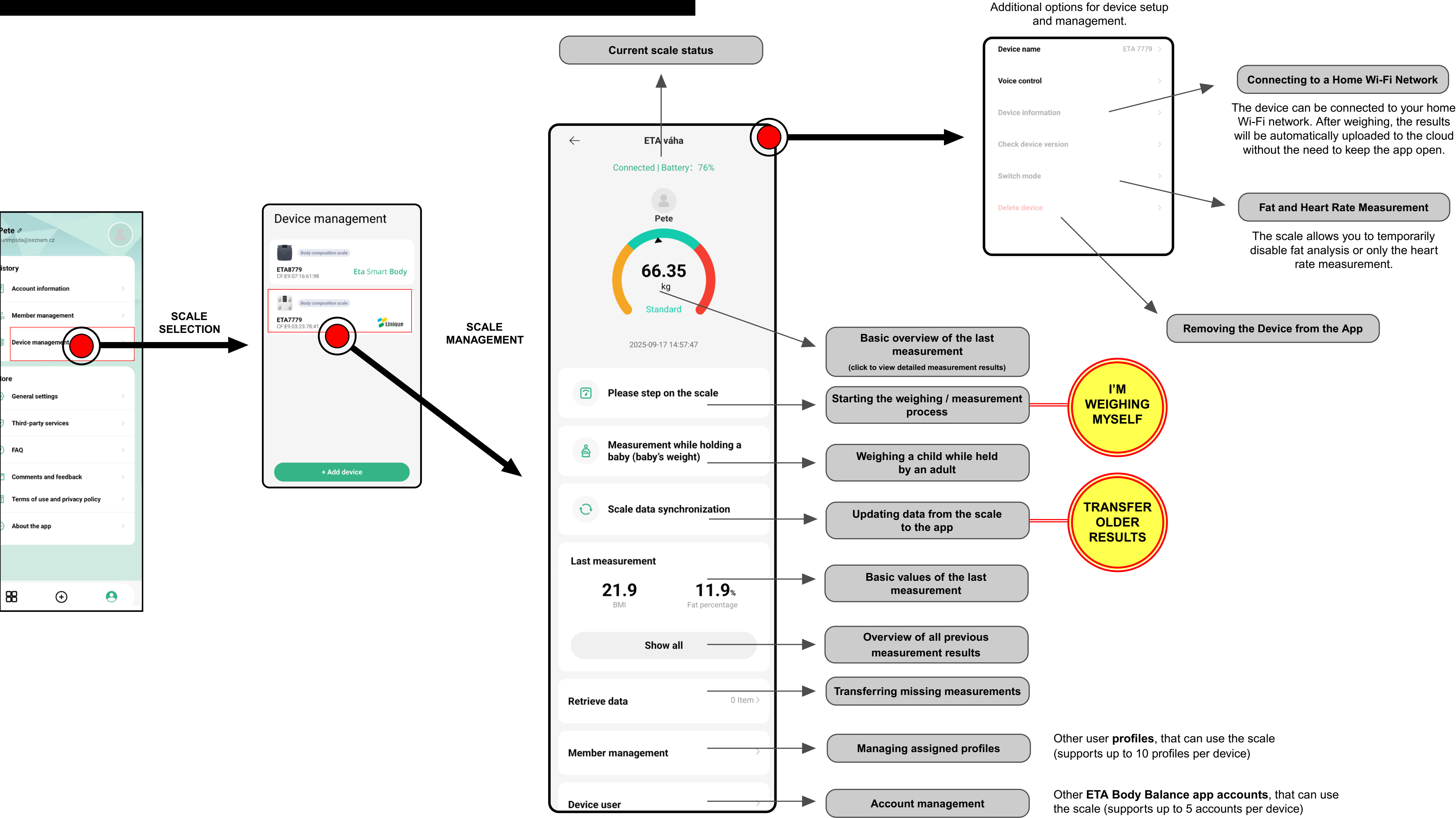USER MANUAL FOR THE Eta Body Balance APPLICATION

Dear Customer,

Thank you for your trust in ETA products. This manual is intended for all owners of the ETA8779, scale who wish to make full use of its potential with the Eta Body Balance application.

For more information, advice, and tips regarding your scale, please visit our website:

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These web pages are regularly updated and contain the latest information about the product.

If you have any questions about the operation of the application or experience any issues, please first visit the above-mentioned website. If the information provided there does not resolve your issue, you can contact us at the following email address: **apps@eta.cz**

1. APPLICATION INSTALLATION

You can find the application in the Apple Store or Google Play application databases, or you can use the following QR code for faster access:

Download and install the application on your mobile phone (or tablet).



COMPATIBILITY	
Mobile device operating system requirements:	Android 7.0 or higher
	IOS 12.0 or higher

1

2

3

2. ACCOUNT REGISTRATION AND LOGIN

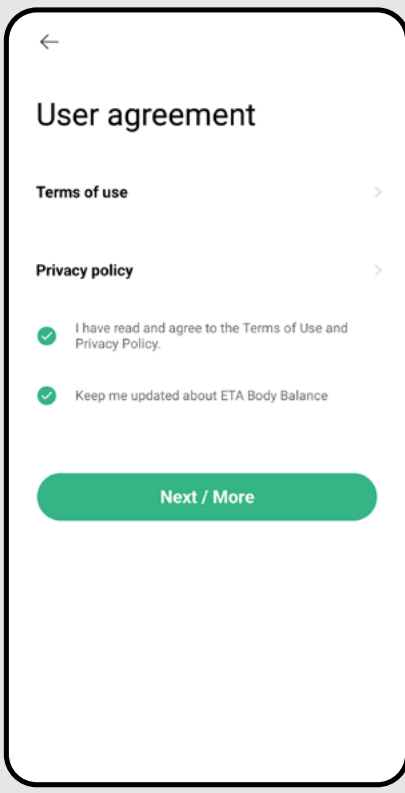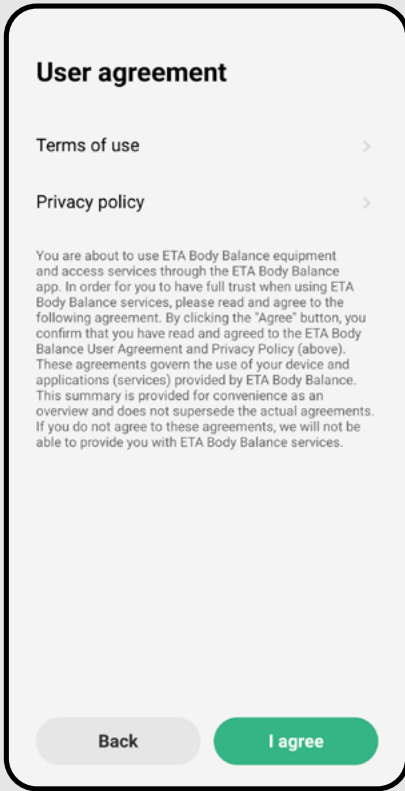- Once the application is downloaded, open it. A legal notice will appear, requiring your consent to the terms of use and privacy policy.
- Click "**I agree**" to proceed to the application.



- Next, the main **login** and new account **registration screen** will appear.
- To log in, click "**Login**" enter your account name (email) and password, then confirm by clicking „**Login**“.
- If you do not have an account yet, continue by selecting „**Registration**“.

- Confirm the necessary information again and proceed by clicking „**Next / More**“.

4

5

6

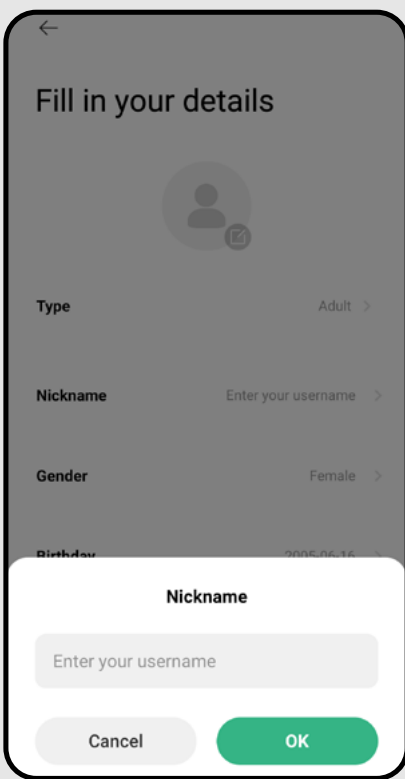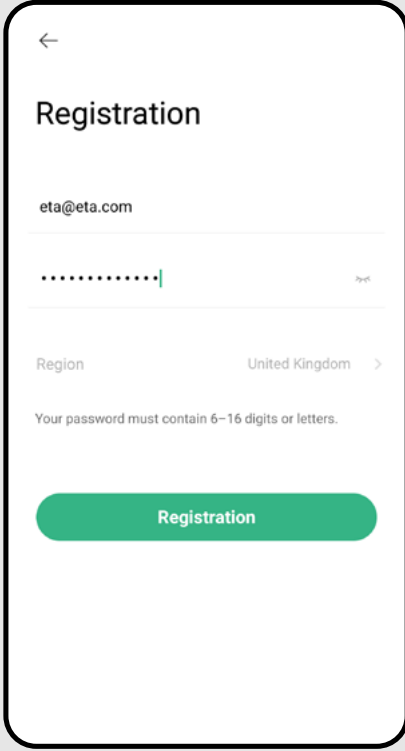
- Enter your email address and create a password for logging into the application.
- Then, select the country in which you will use the device.
- Complete the process by clicking "**Registration**".

- Set the basic parameters of the scale and confirm by clicking "**Next / More**".


- Create a nickname for the main profile and confirm by clicking "**OK**".

7

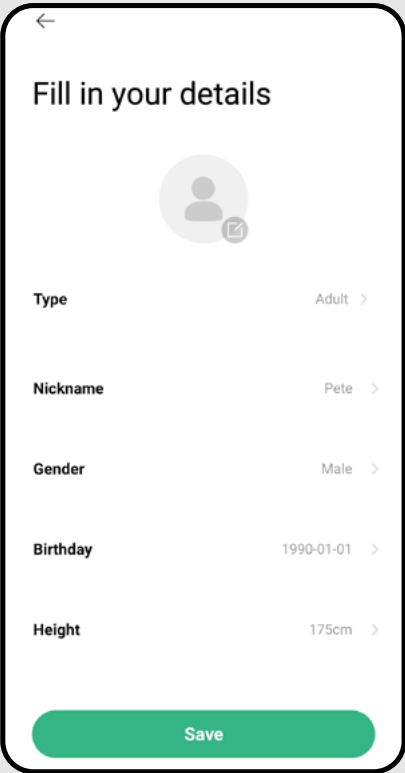- Fill in the main details of your profile, which are necessary for a more accurate calculation of body parameters.
- Confirm by clicking "**Save**".

8

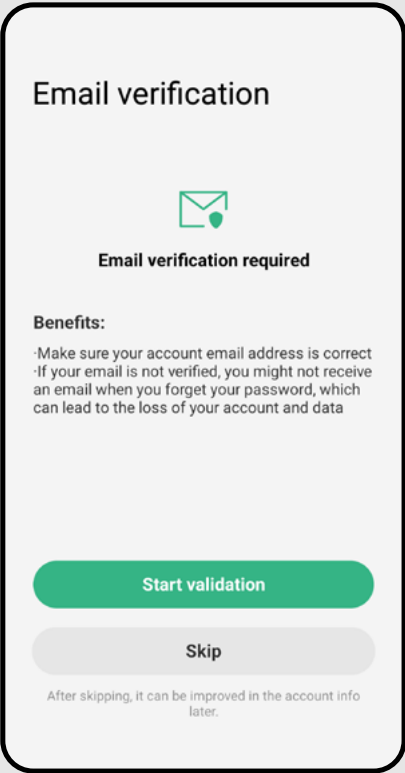
- To enable the full functionality of your account, confirm your email address by clicking "**Start validation**".

9


- 1) Click on the "**Get verification code**" field – a verification code will be sent to your registered email address.
- 2) Enter the verification code into the "**Enter verification code**" field.
- The "**Confirm**", button will then become available to complete the process.

3. DEVICE PAIRING

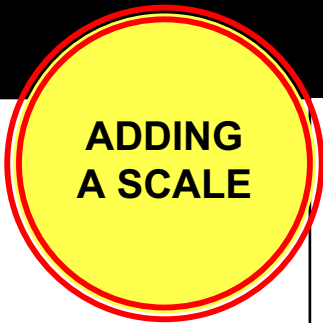
1

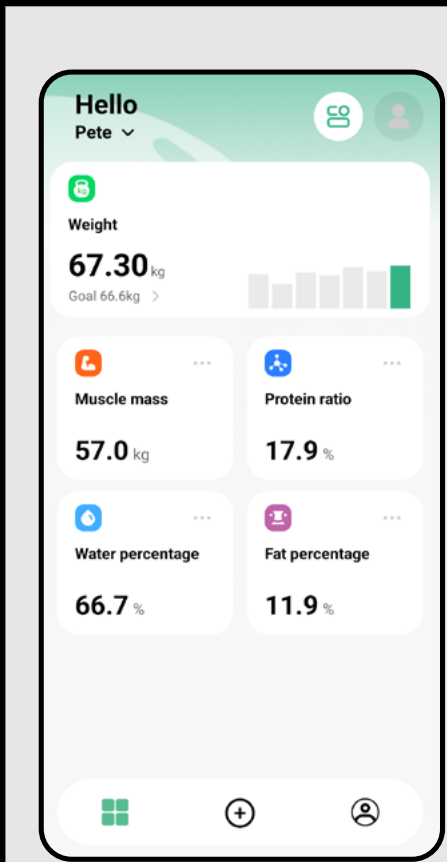
- After successful login, the default home screen will appear.
- To pair the application with the scale, you have several options — see the **DEVICE AND USER MANAGEMENT** page. However, if no device has been paired yet, simply click the "+" icon in the bottom menu.

2

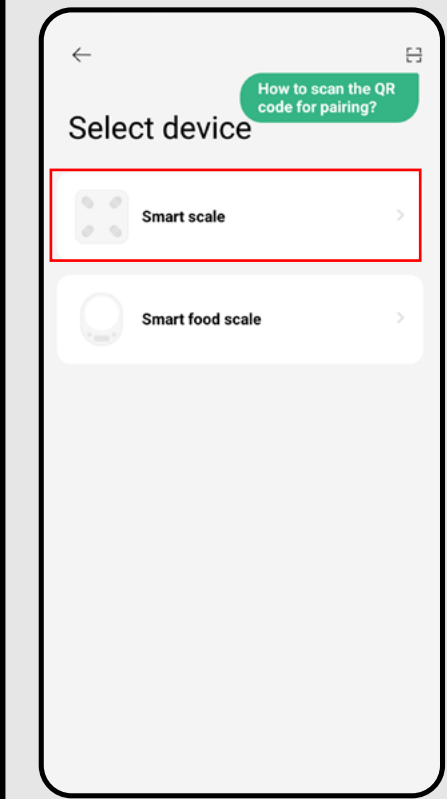
- The application may request permission to access certain features (e.g., Bluetooth). Approve these permissions to ensure proper app functionality.
- Then select the category of scale you want to connect.

3

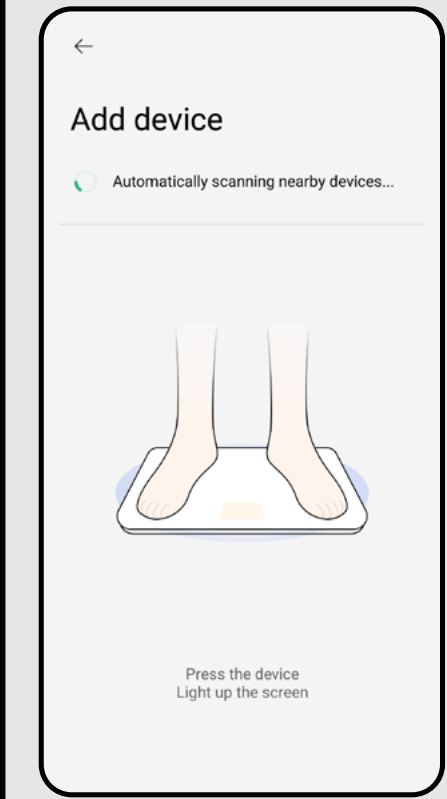
- The app will now automatically start searching for the scale.
- Turn on the scale by stepping on it.

4

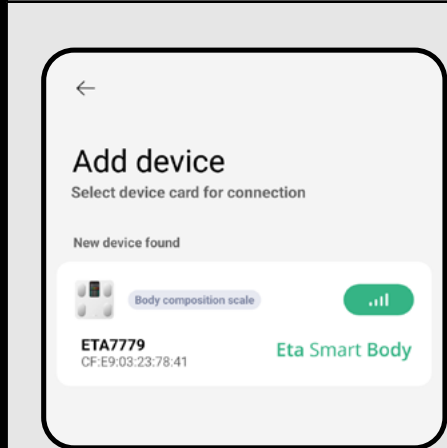
- The application will display all compatible devices found.
- Click on the desired device to start the pairing process.

5

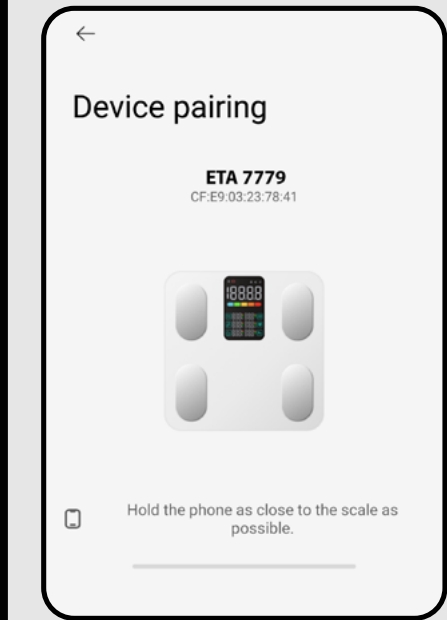
- The pairing process will run automatically.

6

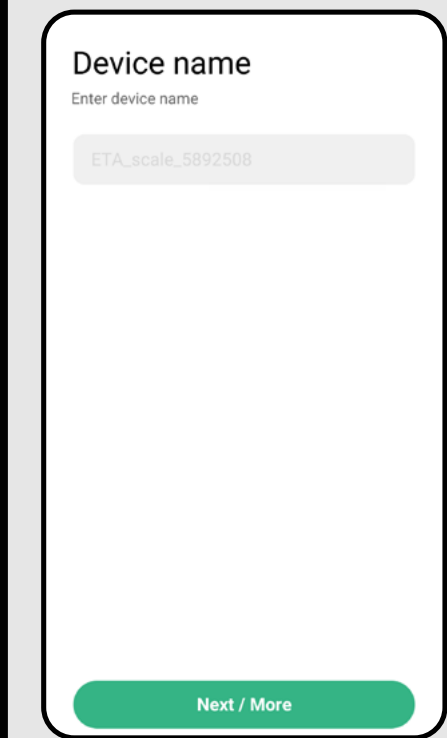- Choose a device name and confirm by clicking "**Next / More**".

7

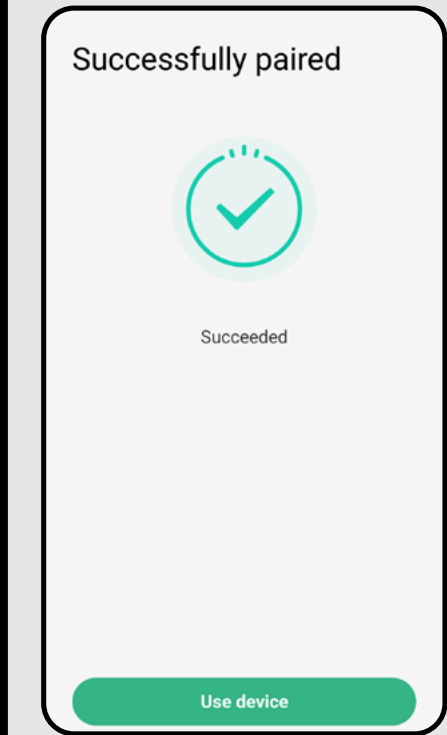- If the pairing is successful, the application will display a notification confirming the connection.
- Finish the process by clicking "**Use device**".

4. BASIC ADVICE AND TIPS FOR USE

Using the ETA scale together with the Eta Body Balance app is intuitive, and during our testing, we did not encounter any significant issues or difficulties you might face during use. However, smart devices in general carry a higher risk of incompatibility with certain mobile devices or potential issues with wireless connectivity. Therefore, we have prepared a set of basic tips and recommendations to help you resolve any issues that may occur when using smart devices.


1. Smart Device (Mobile Phone or Tablet):

- a) Always use the latest version of the app downloaded from official sources.
- b) Make sure your device runs a supported version of Android or iOS.
- c) Ensure that the Bluetooth function is enabled.
- d) Keep your internet connection active so the app can properly communicate with the cloud server.
- e) Grant the app the required permissions, which may vary depending on your phone's operating system. If you experience problems, check whether your phone has any system-level settings that restrict apps from using certain permissions. In such cases, you will need to create an exception for the app.
- f) Communication between the scale and your phone may be disrupted by certain programs (e.g., VPNs, firewalls, etc.) installed on your device. If you experience pairing or operational issues, disable these programs temporarily.
- g) If problems persist, try pairing the scale with another phone to rule out a device-specific issue.

2. Scale Unit

- a) Always make sure you are using sufficiently charged batteries.
- b) Keep the contact surfaces of the scale clean at all times. This minimizes the risk of incorrect measurements. Step onto the scale barefoot and slightly moisten your feet to ensure better conductivity.


- Try creating a new user account to eliminate issues related to account activation.
- During pairing and operation, keep your smart device close to the scale (ideally within approximately 2 meters).

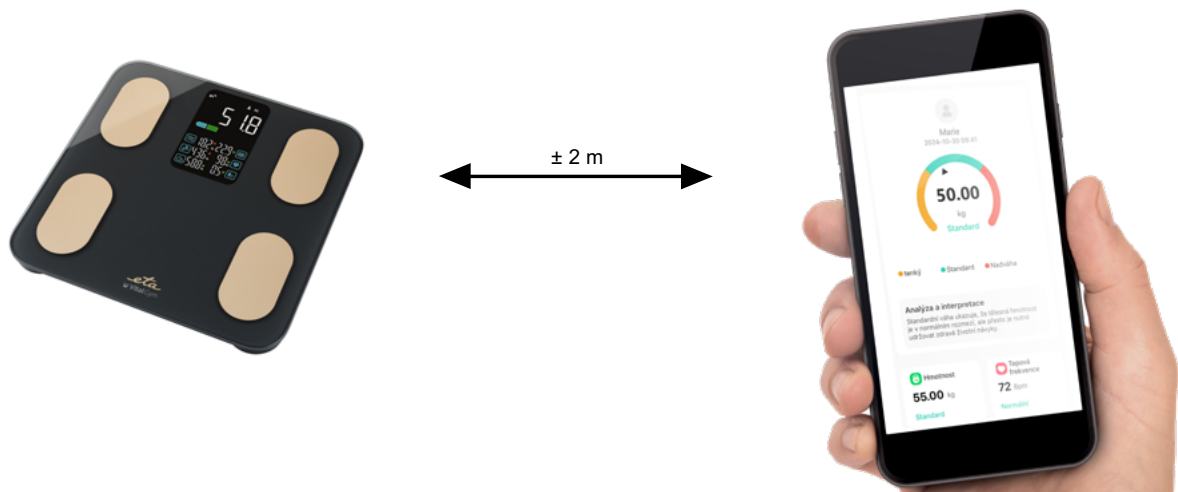5. WEIGHING AND ANALYSIS

WEIGHING

1

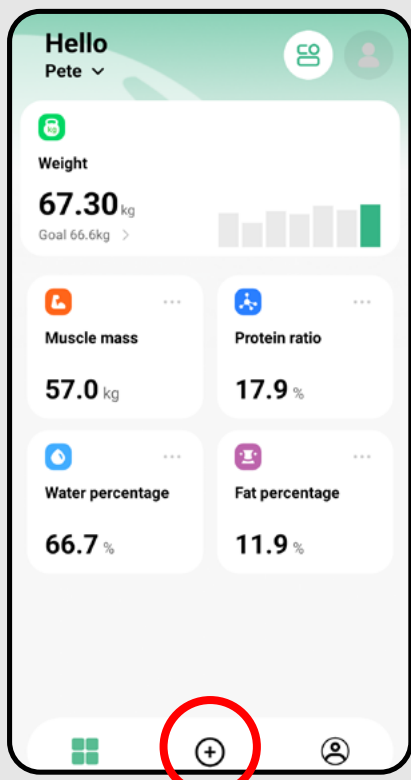
- To start the measurement process, click the "+" icon in the bottom bar.

Note:


If you only want to transfer already measured data, go to Data Synchronization (see section **SCALE SETTINGS AND MANAGEMENT**).

2

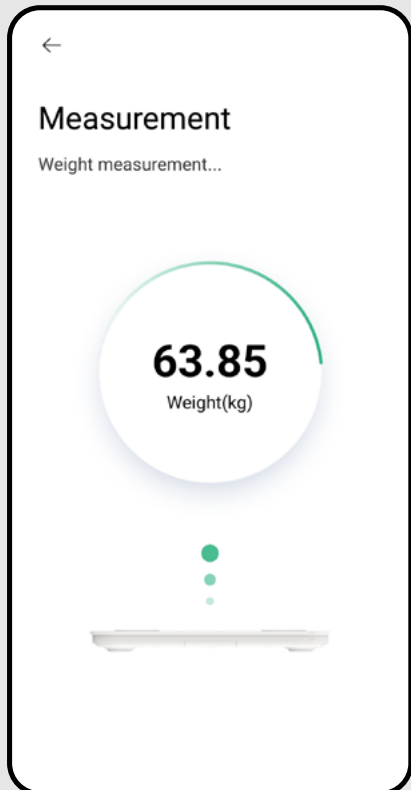
- The first measurement phase – weight measurement – will begin.

3

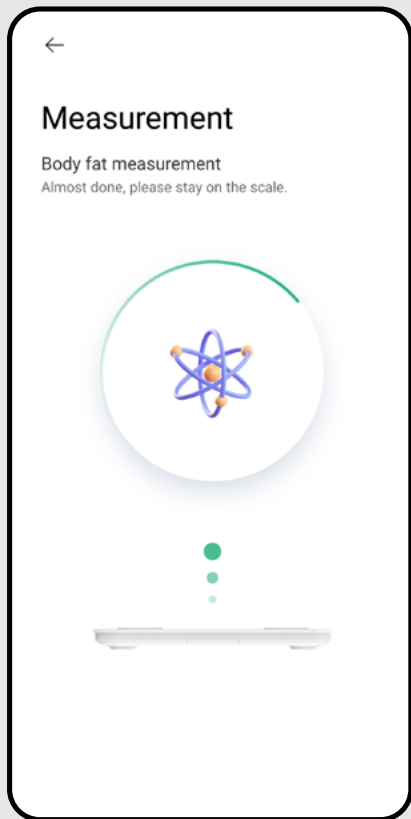- Then the scale will start analyzing body parameters.


4

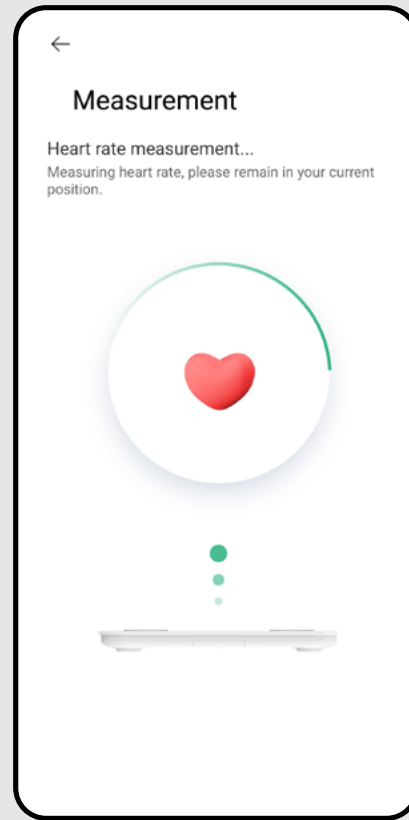- In the final phase, it will measure heart rate.

6. MEASUREMENT RESULTS

1

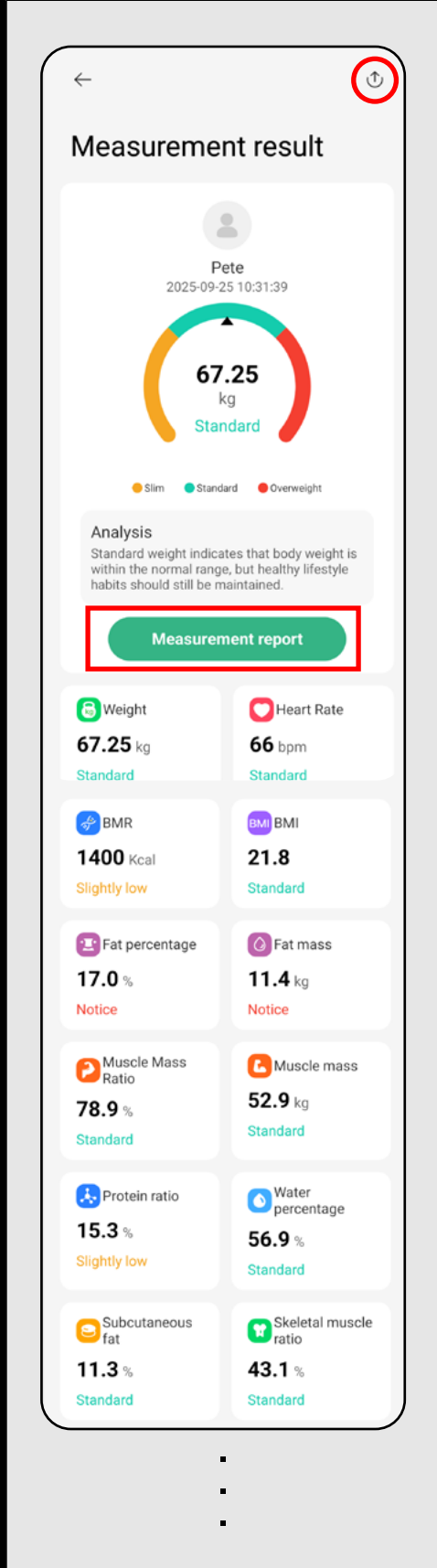• After successful analysis, a screen will appear showing the results of all measured parameters.

• For a more detailed analysis of the measurement results, click **"Measurement report"**.

Note:

• If the analysis does not complete correctly, the application will display only weight and BMI.

SHARING

• To share the measurement results as an image (JPG), click the share icon in the upper right corner.


2

• The displayed measurement report provides a more detailed view of the results.

SUMMARY RESULTS OVER A LONGER PERIOD

- The home screen displays the results of the most recent measurement. By clicking on the selected parameter on this screen, you can view summary results over a longer period.
- The results in the charts can be filtered by individual parameters.

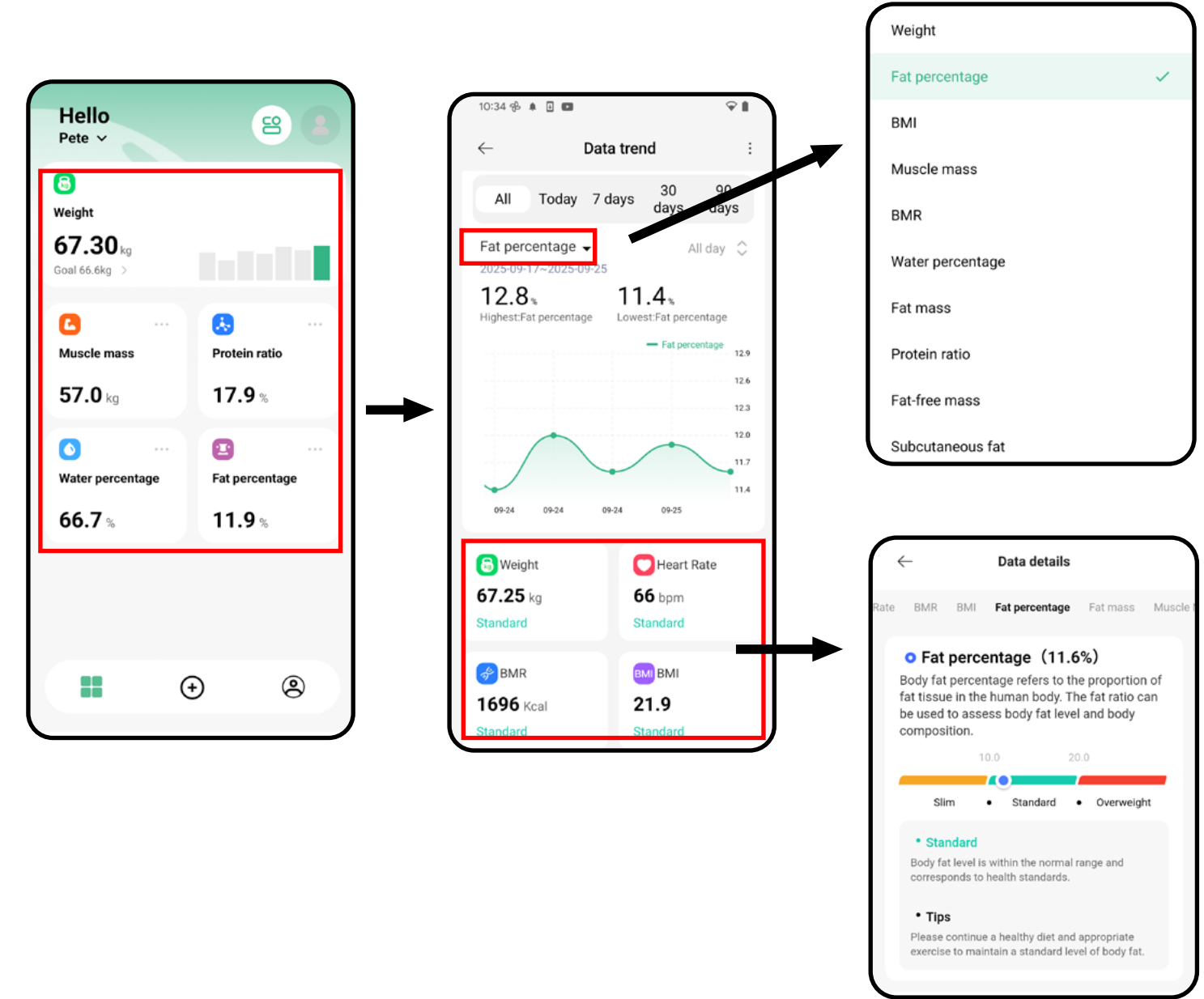
- When you open a specific parameter, brief information and a graphical representation of the measurement result will be displayed.


NOTES

- Please note that the measured values are only approximate and may differ from the results of actual medical analyses. Since measuring the human body is very complex and costly, different manufacturers use various indicative measurement methods. As a result, even professional devices may produce differing results. For this reason, it is recommended to use a single device and monitor the trend of each parameter over time.
- If the scale connection fails, make sure Bluetooth is enabled on your phone. Alternatively, refer to the instructions in Chapter 4.
- If the measurement or analysis fails, try repeating the weighing/measurement process. Alternatively, refer to the instructions in Chapter 4.

From the **Summary Results**, section, where you can find **Data Trends**, you can access the list of performed measurements. Click on  and continue as described below.

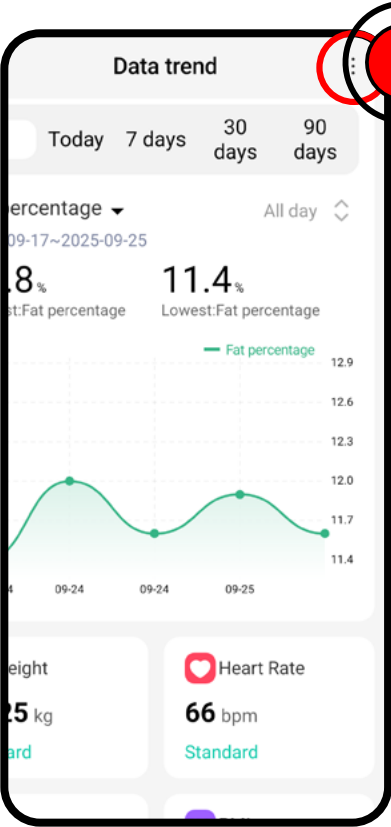LIST OF MEASUREMENT RECORDS

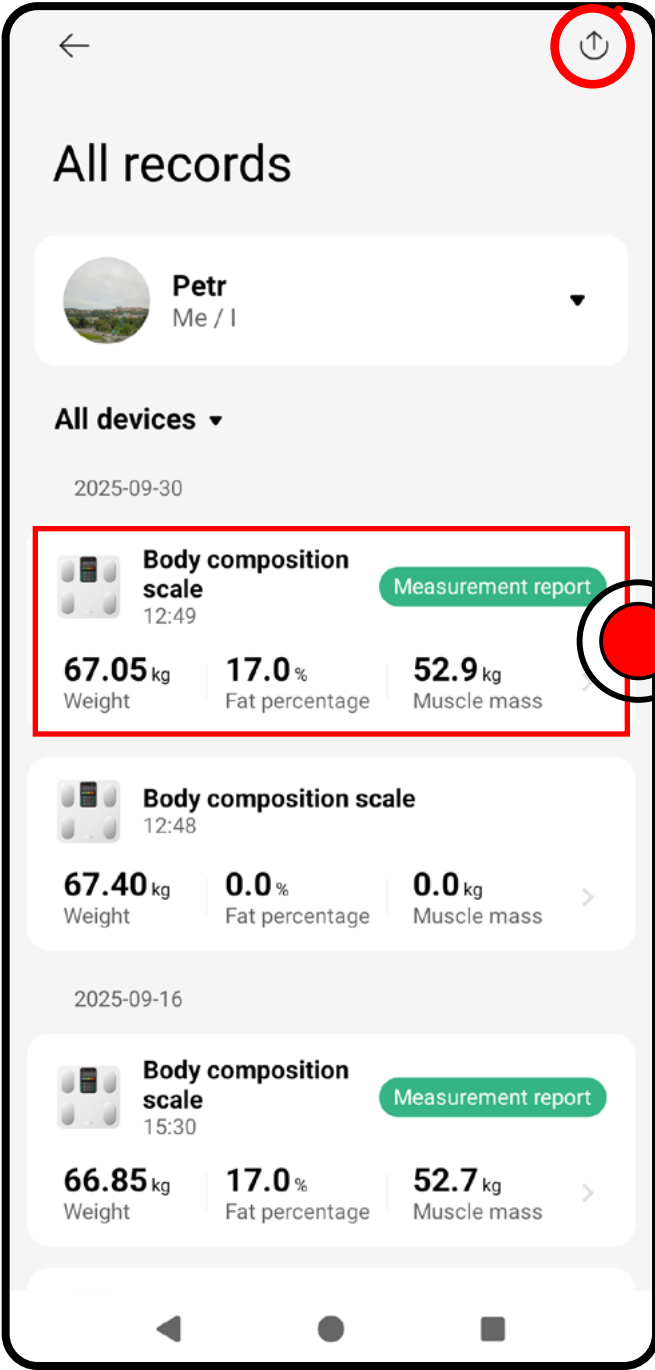Share a basic overview of the list of measurement items in the form of an Excel file.


Share the result of an individual measurement in the form of a JPG image file.

Share the detailed measurement result in the form of a URL link.



BASIC RESULTS OF A SINGLE MEASUREMENT



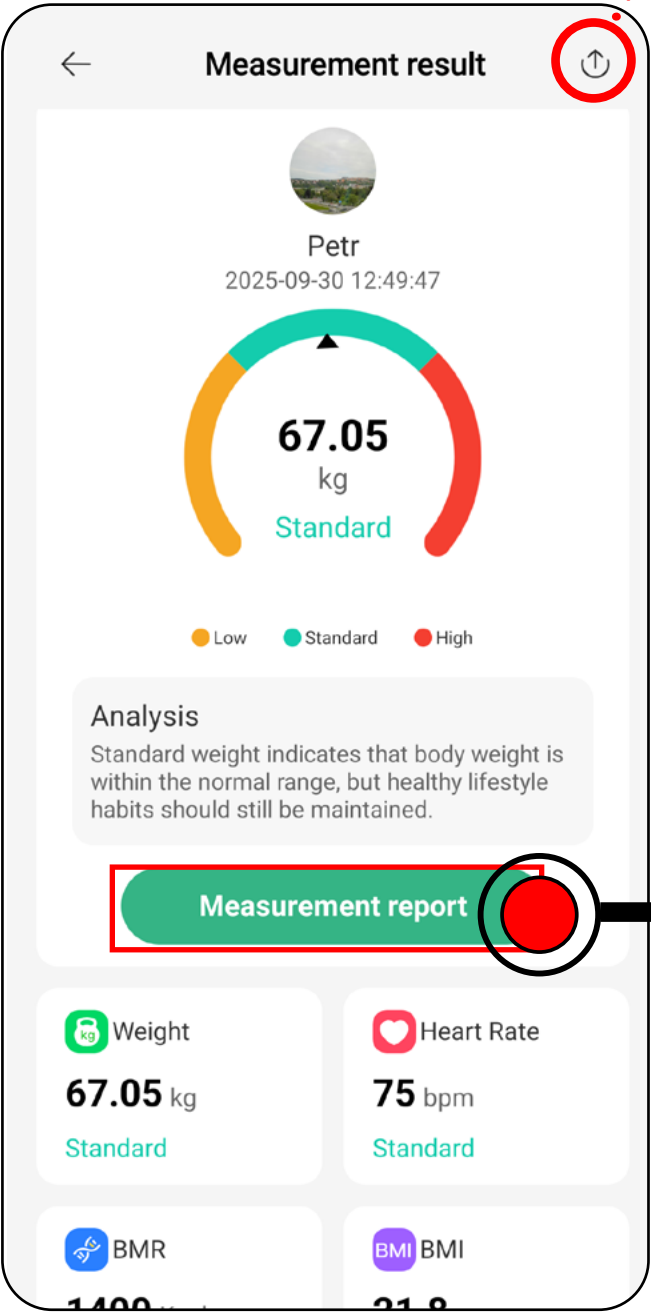
DETAILED RESULTS



8. THIRD-PARTY APPLICATIONS

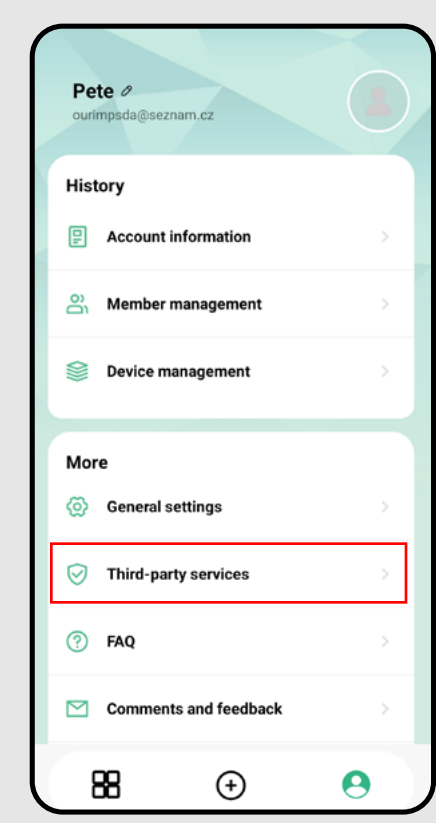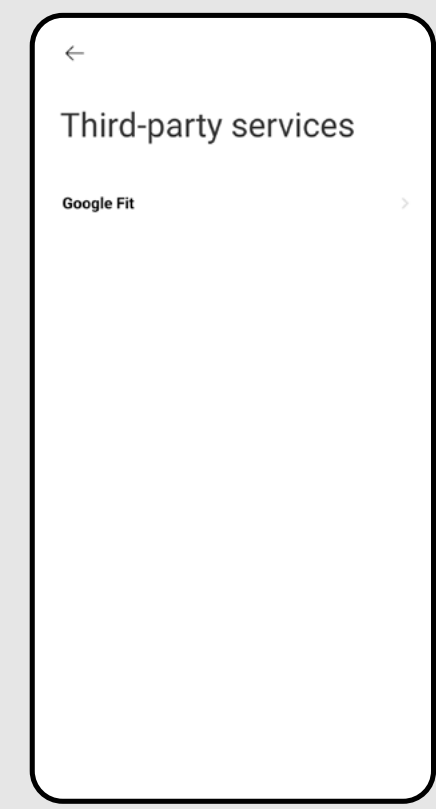1		• The application allows sharing selected result parameters with third-party applications. • After connecting, the respective application will download the measurement results of the selected parameters. • Click "Third-party services".
2		• For Google, this refers to the Google Fit application. • For Apple, this refers to the Apple Health application. • Click the respective program and follow the on-screen instructions. Note • The application will display only the option supported by your phone's operating system.

9. QUESTIONS AND ANSWERS

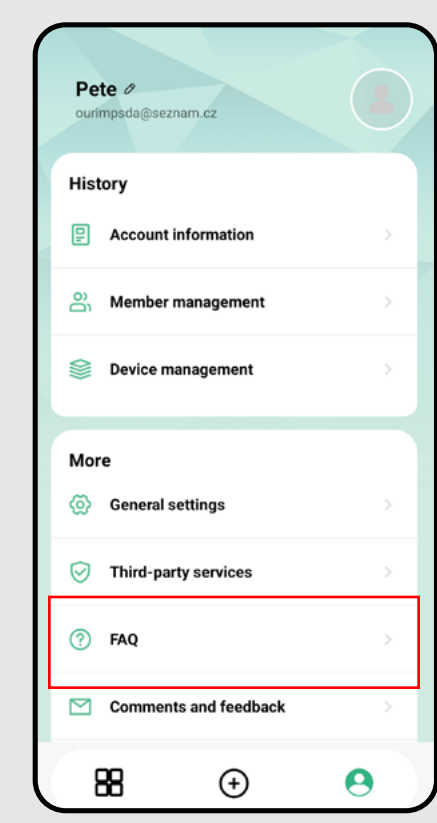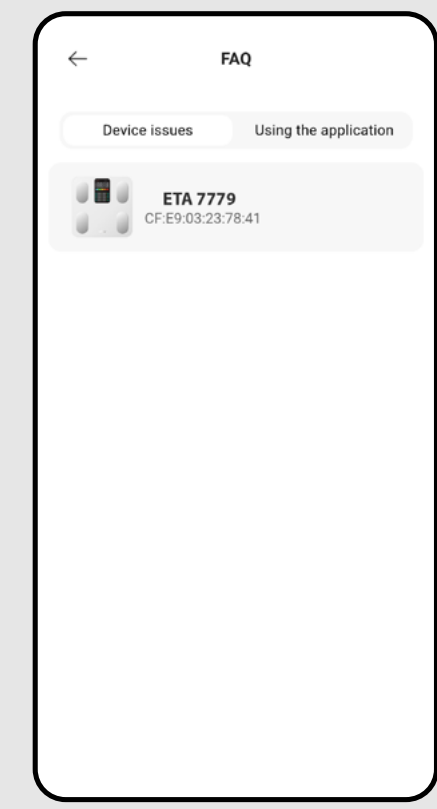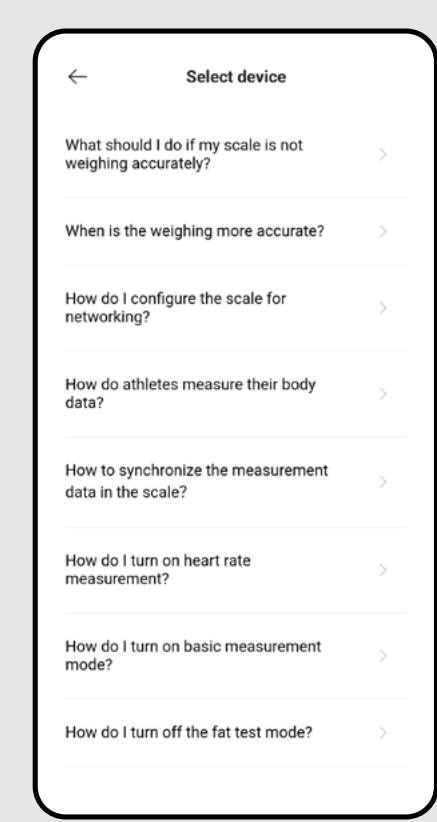1		• The application includes a section with basic questions and answers. • Click "FAQ".
2		• If you have paired multiple devices, all of them will be displayed on this screen. • Click to select a specific model.
3		• Select the question you are interested in.

10. FEEDBACK

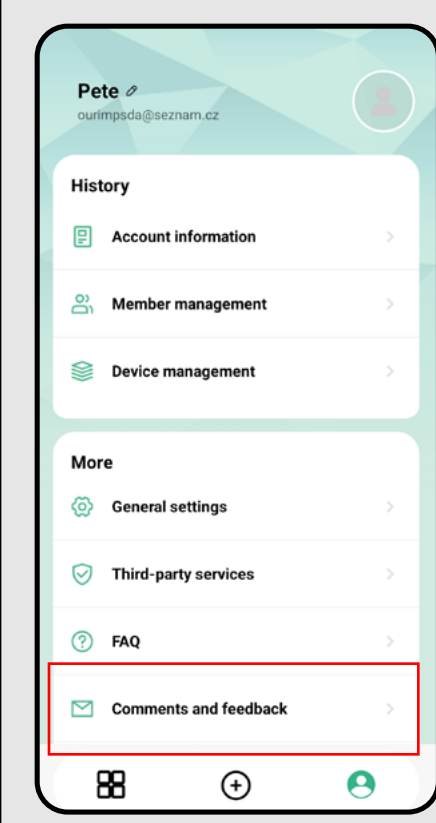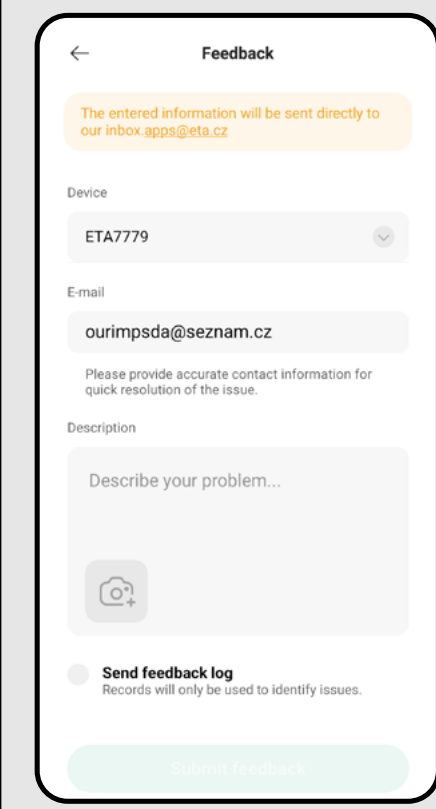1		• If you need assistance or would like to provide feedback, click "Comments and feedback".
2		• Describe your issue or inquiry and send the message by clicking "Submit feedback".

In case of application issues, you can also contact us directly at **apps@eta.cz**

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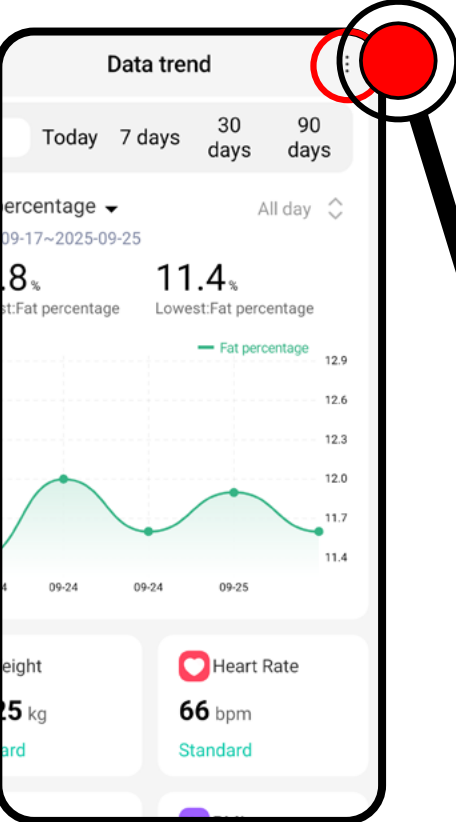
2025-09-30

Body composition scale 12:49

67.05 kg Weight 17.0% Fat percentage 52.9 kg Muscle mass

2025-09-16

Body composition scale 15:30

66.85 kg Weight 17.0% Fat percentage 52.7 kg Muscle mass

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67.05 kg Standard

Low Standard High

Analysis
Standard weight indicates that body weight is within the normal range, but healthy lifestyle habits should still be maintained.

Measurement report

Weight 67.05 kg Standard Heart Rate 75 bpm Standard BMR BMI

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← In-depth analysis ⌵

ETA Body Balance Body composition analysis report

Physical testing Name: Petr Gender: Male Date: 30-09-2025 12:47:47 Age: 39 Height: 175cm

Body composition

Weight	67.0	67.0	67.0
Height	175	175	175
Body mass index	21.4	21.4	21.4
Weight Standard	Water percentage Standard	For men	For men

Comprehensive physical condition

BMI	21.4	21.4	21.4
Heart Rate (bpm)	75	75	75
BMR	1400	1400	1400
BMI	21.4	21.4	21.4

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